


AGILECATH™

Epidural Catheter for Pain Management

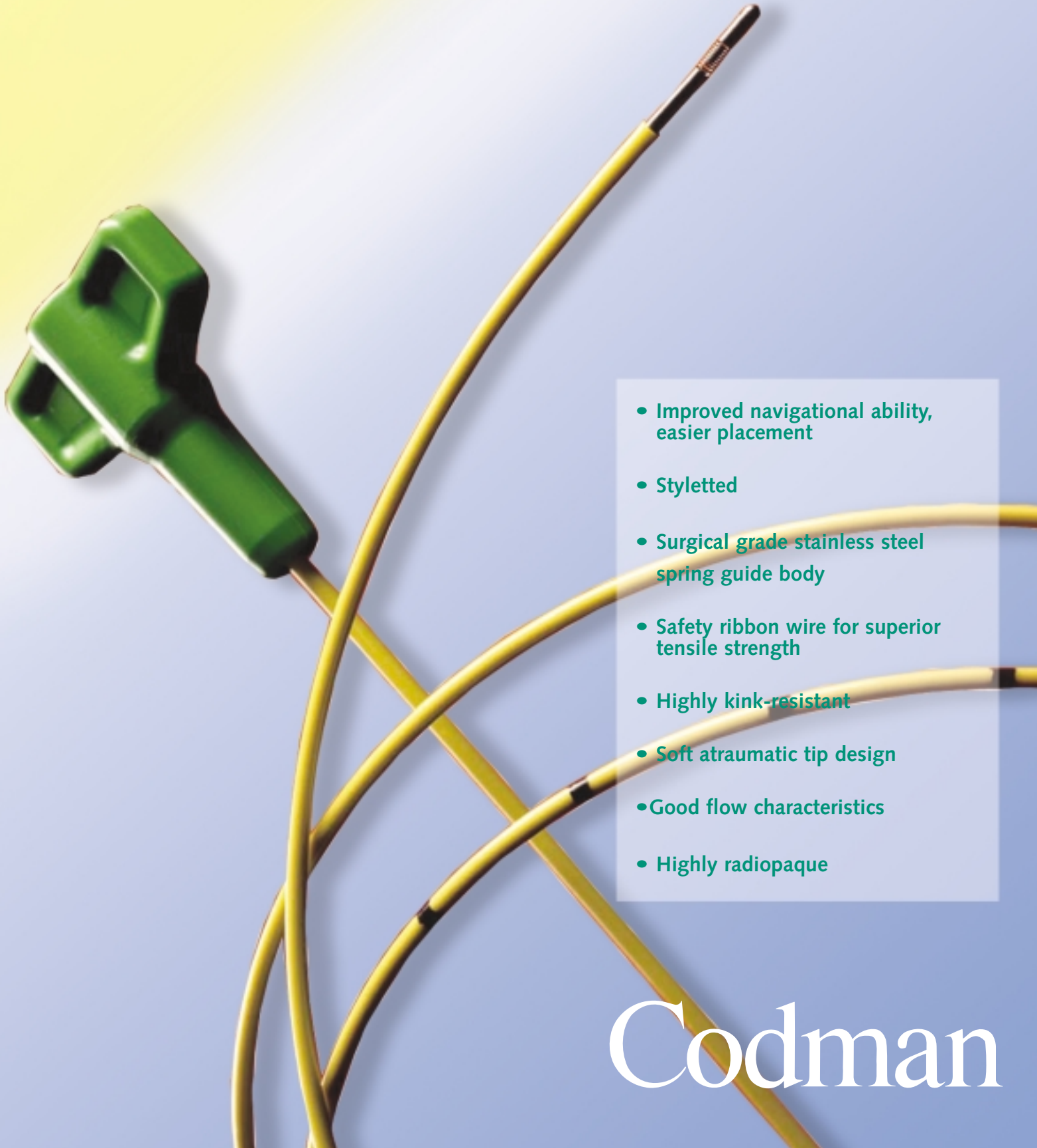

- Improved navigational ability, easier placement
- Styletted
- Surgical grade stainless steel spring guide body
- Safety ribbon wire for superior tensile strength
- Highly kink-resistant
- Soft atraumatic tip design
- Good flow characteristics
- Highly radiopaque

Codman

AGILECATH™ Epidural Catheter for Pain Management

Epidural catheter placement can be challenging. The AGILECATH addresses problems inherent to catheter placement—column-strength and steerability—to make catheter placement easier and more precise.

- A smooth, durable, polymer exterior and spiral-wound stainless steel helix provide high compression strength to resist collapsing and promote positive patency over extended use.
- Improved column-strength stylet allows easier navigation through the epidural space.
- An epidural introducer helps make insertion smoother and reduces the potential for catheter shearing.



AP-01900 AGILECATH Epidural Catheter Introducer

- A stainless steel helical wire identical to the distal soft-tip spring-wire guide allows solutions injected through the internal lumen to be easily dispersed through the expanded windings at the uncovered distal tip.
- Greatly enhanced catheter radiopacity.
- A simple-to-use snaplock catheter/syringe adapter can be secured positively to the catheter without collapse of the lumen.



Codman
a Johnson & Johnson company

For more information, contact your Codman Drug Delivery Representative

325 PARAMOUNT DRIVE • RAYNHAM, MA 02767 • (508) 880-8100
FOR PRODUCT INFORMATION, CALL: (800) 225-0460 • www.codmanjnj.com

Ordering Information

AGILECATH Epidural Catheter Systems

Product No.	Description	Quantity / Case
AP-01900	AGILECATH Epidural Catheter Introducer	10
AP-01931	One: 31cm length 19 Ga. AGILECATH Epidural Catheter with Stylet, snaplock Adapter One: Epidural Catheter Introducer: 17 Ga. Radiopaque Polyurethane with Integral Suture Wing One: Epidural Needle: 17 Ga. x 3-7/8" (9.84cm) TW with cm markings	10
AP-01960	One: 60cm length 19 Ga. AGILECATH Epidural Catheter with Stylet, snaplock Adapter One: Epidural Catheter Introducer: 17 Ga. Radiopaque Polyurethane with Integral Suture Wing One: Epidural Needle: 17 Ga. x 3-7/8" (9.84cm) TW with cm markings	10
AP-01990	One: 90cm length 19 Ga. AGILECATH Epidural Catheter with Stylet, snaplock Adapter One: Epidural Catheter Introducer: 17 Ga. Radiopaque Polyurethane with Integral Suture Wing One: Epidural Needle: 17 Ga. x 3-7/8" (9.84cm) TW with cm markings	10
AP-01931-PK	Full Procedure Kit, 31cm length AGILECATH Epidural Catheter	10
AP-01960-PK	Full Procedure Kit, 60cm length AGILECATH Epidural Catheter	10
AP-01990-PK	Full Procedure Kit, 90cm length AGILECATH Epidural Catheter	10