
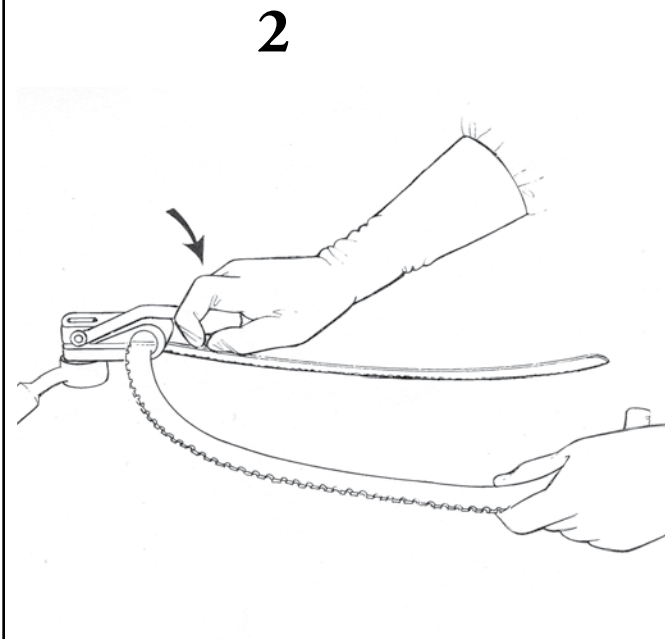


Set-Up

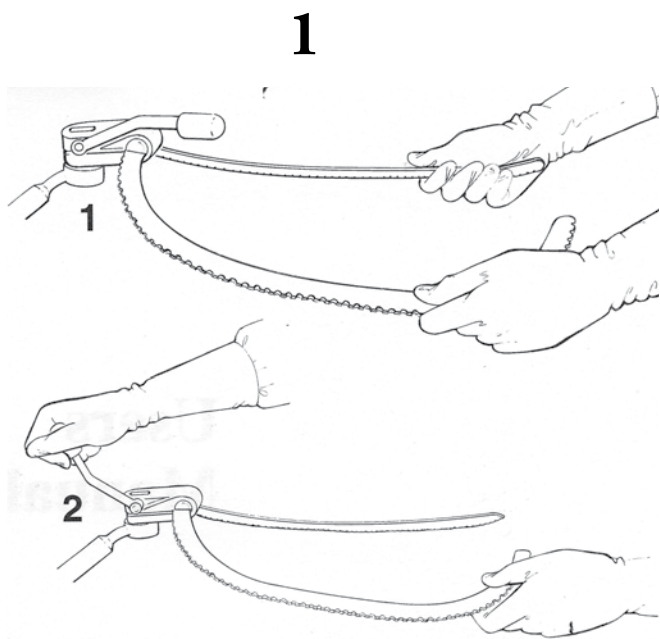

Insert the BOOKWALTER WISHBOOK support arm into the post coupling and secure the frame over the patient (1). Adjust the frame arms into the desired position. Lock the arms in position by pushing the securing handle back (2).

To re-adjust the frame, unlock the securing handle.

Ordering Information

- 50-4648 - BOOKWALTER WISHBOOK Adjustable Frame
- 50-4649 - BOOKWALTER WISHBOOK Adjustable Frame - Small




- NECESSARY FOR COMPLETE SETUP:**
- 50-4581 Bookwalter TablePost
  - 50-4554 Bookwalter Post Coupling
- OPTIONAL:**
- 50-4589 Bookwalter ROTILT™ Ratchet Mechanism


BOOKWALTER™ is a registered trademark of John Bookwalter, M.D.  
 WISHBOOK™ is a trademark of Minnesota Scientific, Inc.  
**TO ORDER CALL (800) 255-2500**



For more information, contact your Codman Sales Representative.  
 325 PARAMOUNT DRIVE • RAYNHAM, MA 02767 • (508) 880-8100  
 FOR PRODUCT INFORMATION, CALL: (800) 225-0460

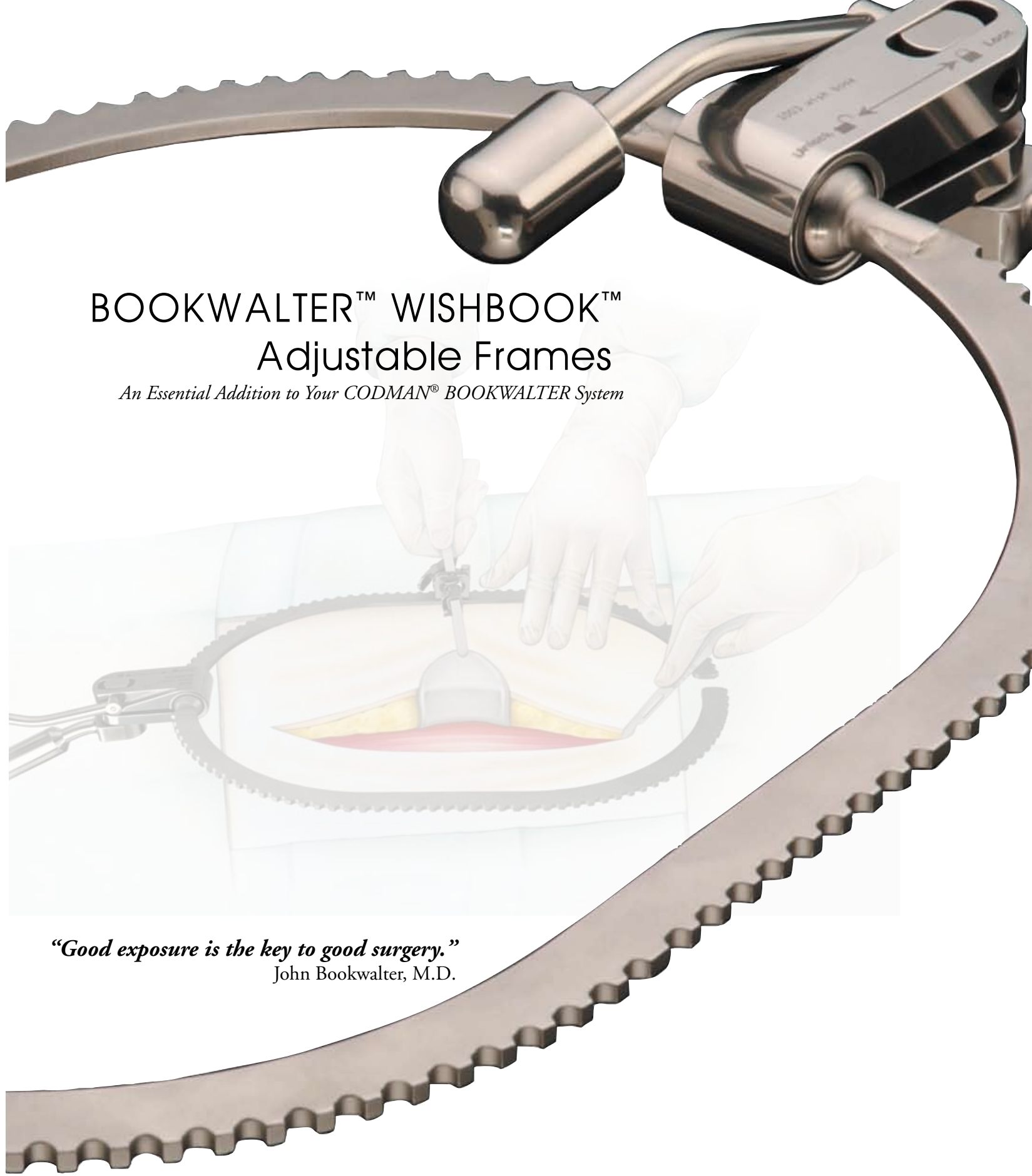

BOOKWALTER™ WISHBOOK™  
 Adjustable Frames  
 An Essential Addition to Your CODMAN® BOOKWALTER System

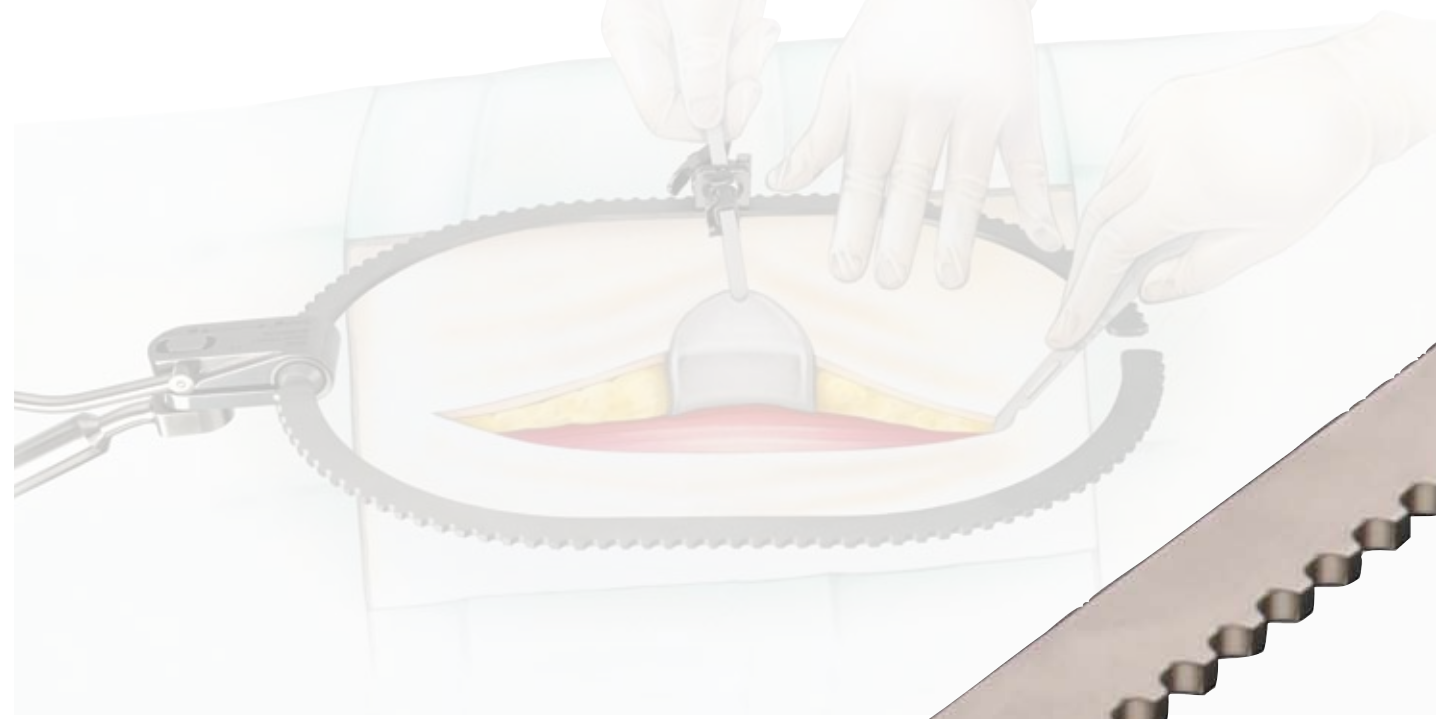

*“Good exposure is the key to good surgery.”*  
 John Bookwalter, M.D.




# BOOKWALTER™ WISHBOOK™ Adjustable Frames

Improve the performance of your CODMAN BOOKWALTER System.

The New BOOKWALTER WISHBOOK Adjustable Frames reduce critical set-up time during emergency, trauma and various other open surgical procedures.

The BOOKWALTER WISHBOOK is designed to improve the function and flexibility of your Bookwalter Retractor System. With the BOOKWALTER WISHBOOK, frame arms can be placed exactly where you want them and can easily be adjusted anytime during the procedure.

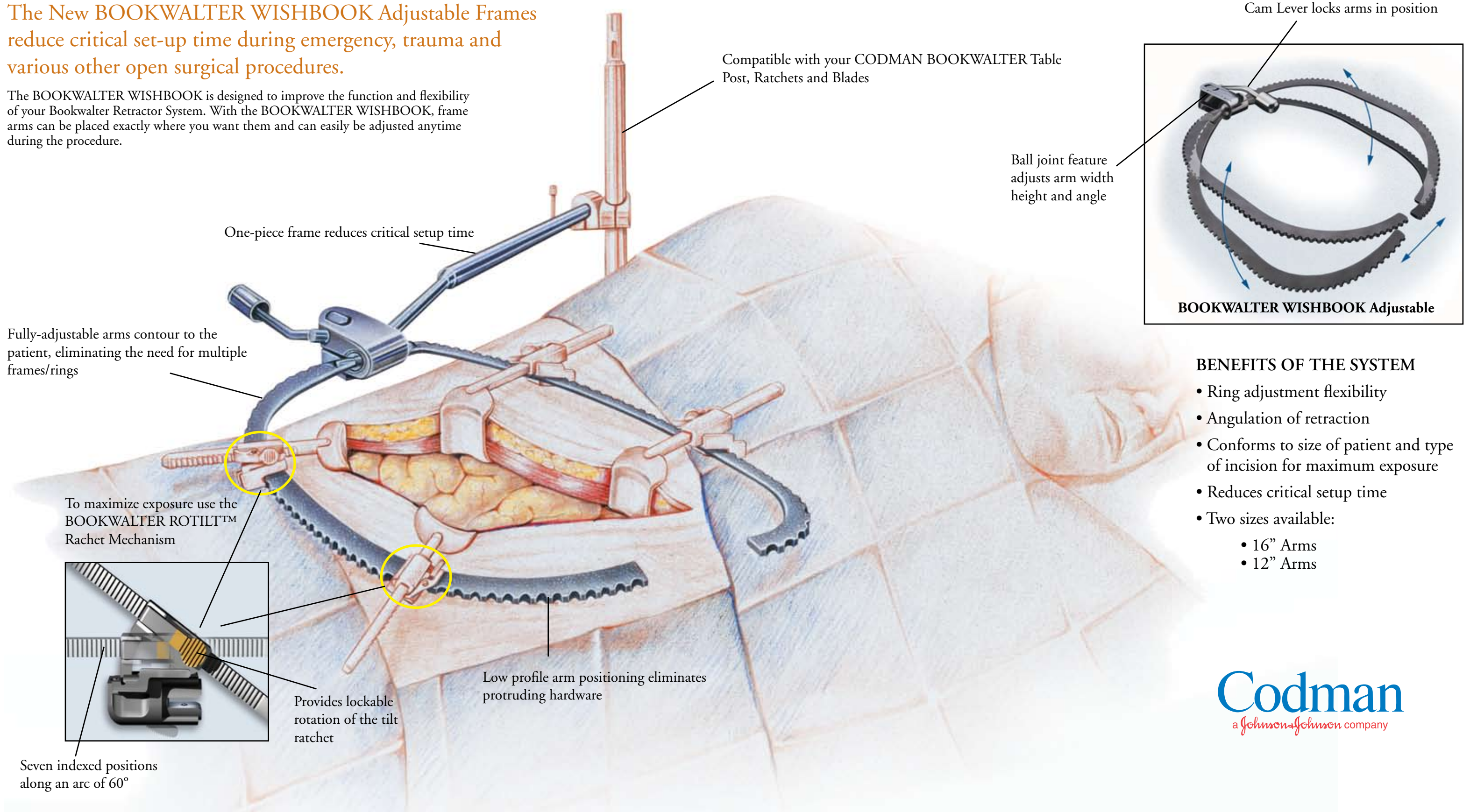

## BENEFITS OF THE SYSTEM

- Ring adjustment flexibility
- Angulation of retraction
- Conforms to size of patient and type of incision for maximum exposure
- Reduces critical setup time
- Two sizes available:
  - 16" Arms
  - 12" Arms

**Codman**  
a Johnson & Johnson company